

## CONNECT YOUR DEVICE TO MOTICAM T2 PORTABLE HOTSPOT

Step 1. Set up the Portable Hotspot on your Moticam T2.

Turn on the Motic Hub (fig. 1).

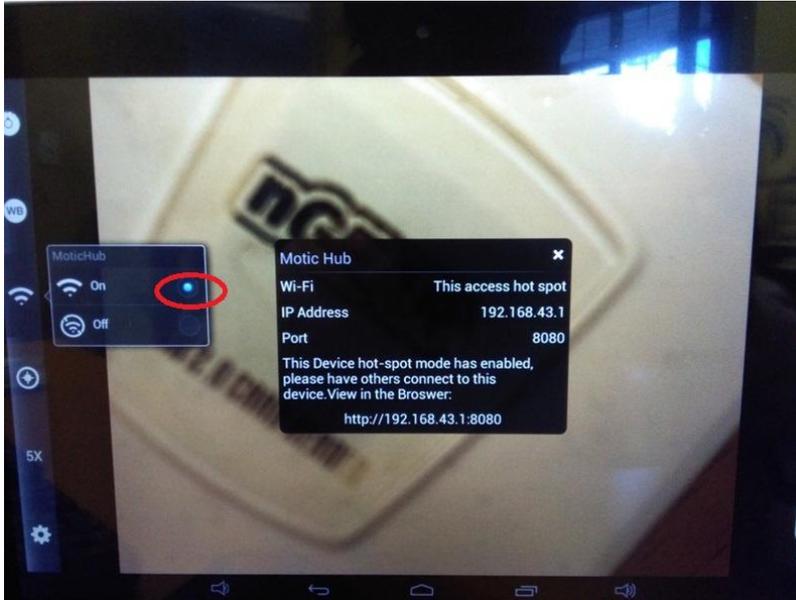


fig. 1

Step 2. You can use either a web browser, MotiConnect app or Motic Image Plus 2.0 ML software to connect to MoticamT2 Portable Hotspot

- Using a web browser:

Type the Moticam T2 IP address and port number (fig. 2) into a web browser address bar (fig. 3). You can get the IP address and port number by clicking on Settings --> IP Address.



fig. 2



fig. 3

- Using MotiConnect app

Open the app, click on Live Camera.

Click on Options --> Add, type in the MotiCam T2 Wi-Fi Network name, IP address and Port number (fig. 4) to your device (fig. 5)



fig. 4

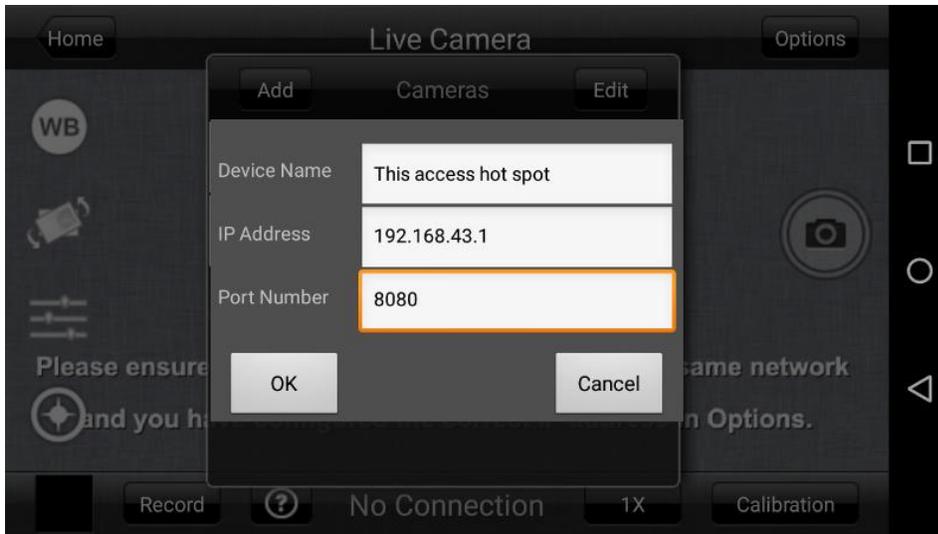


fig.5

Click on OK. Then make sure the camera is chosen (fig. 6)



fig.6

- Using Motic Image Plus 2.0 ML software  
Open the software, click on Capture --> Open, make sure the **Video Device** is Moticam Xtype in the Moticam T2 IP address (fig.7).

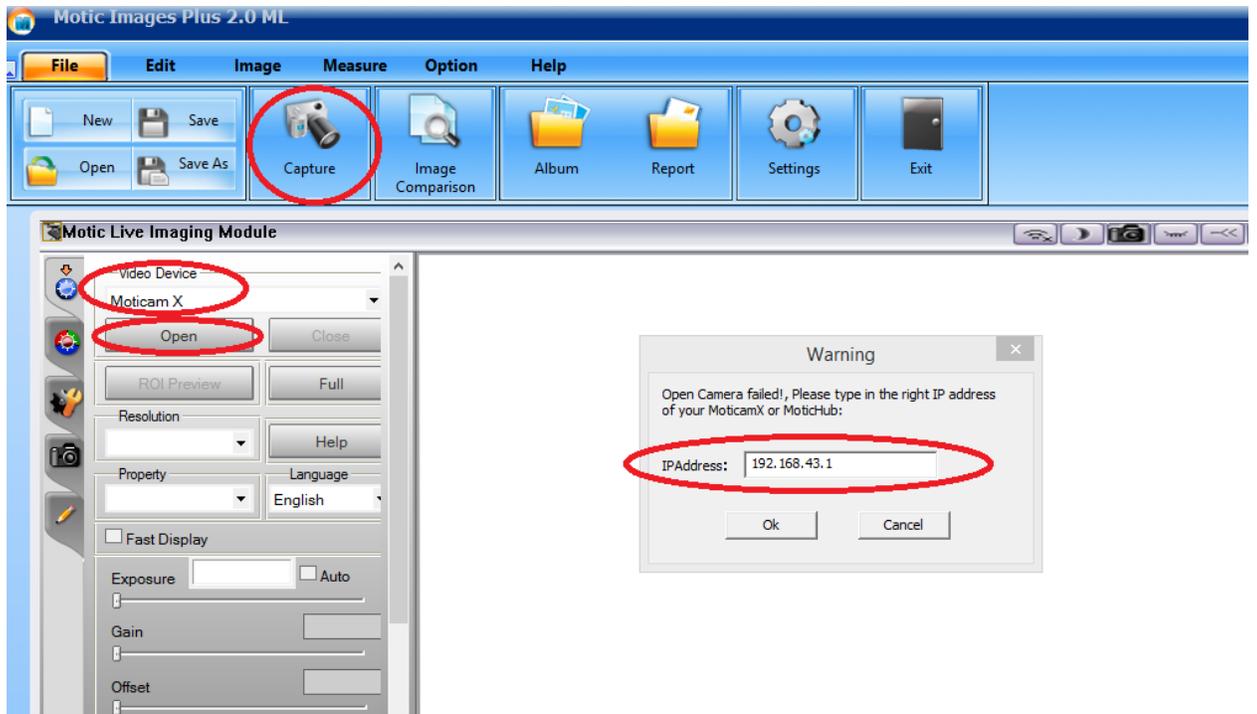


fig. 7